


# Optima Millennium

## *Bucky stand*

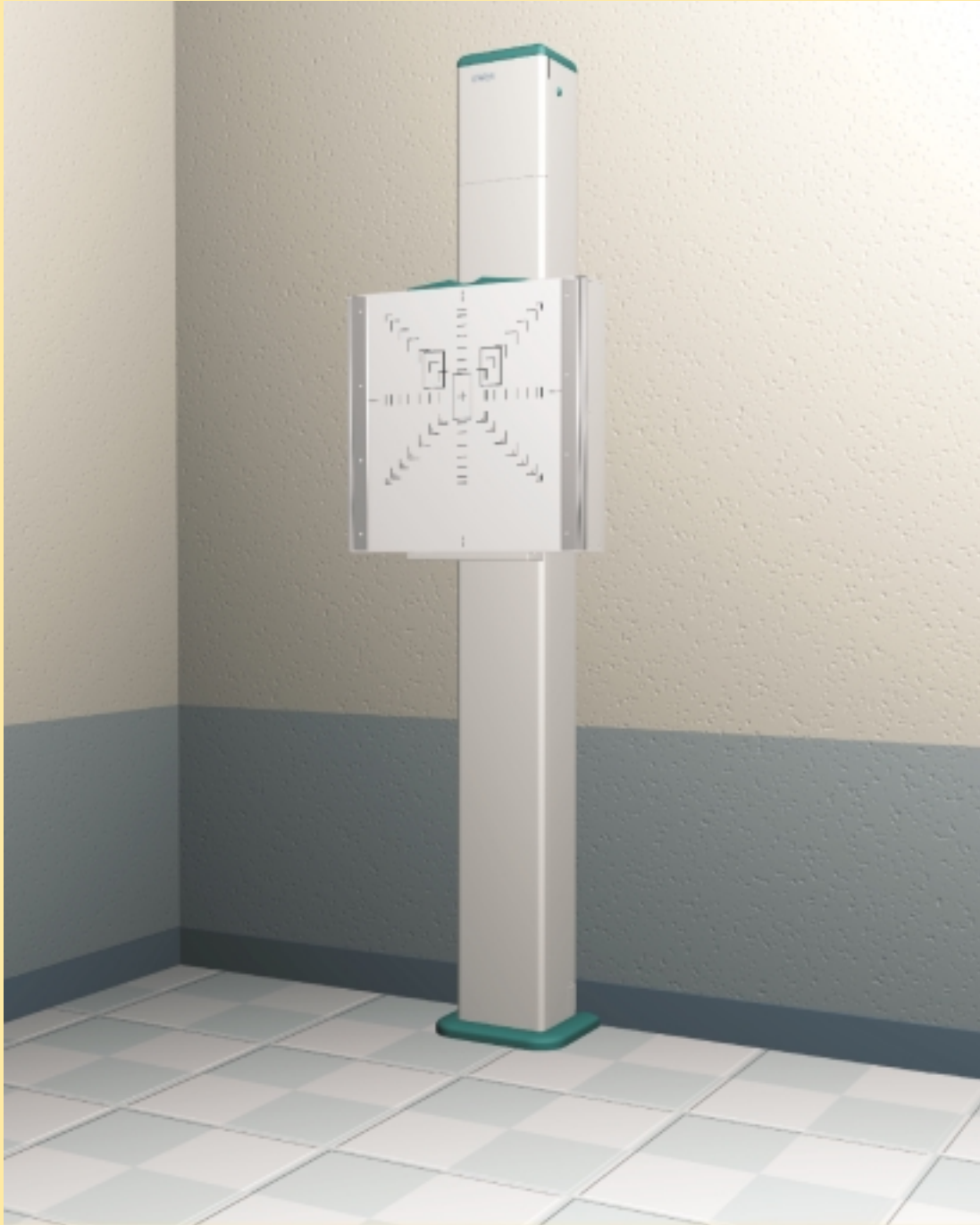



# Optima Millennium

## Wall Bucky

### Description

Wall fixed stand with vertical sliding bucky carriage. Manual moving, weight counterbalancing, and vertical positioning by means of electromagnetics breaks. The bucky carriage can accept a standard bucky mechanism with 35x43 cm cassette tray and grid. The bucky front panel is made of a special thermo stable resin with very low absorption level, a chin-support on top side, and the lateral aluminium rails fitted with grooves for standard accessories. The cassette loading can be done on the right or on the left side. As option it is possible to mount the automatic cassette loading bucky mechanism and the cassette size sensing device for use with an automatic collimation system. All buckys accept a lon chamber for use with an automatic exposure control.

## DIMENSIONS

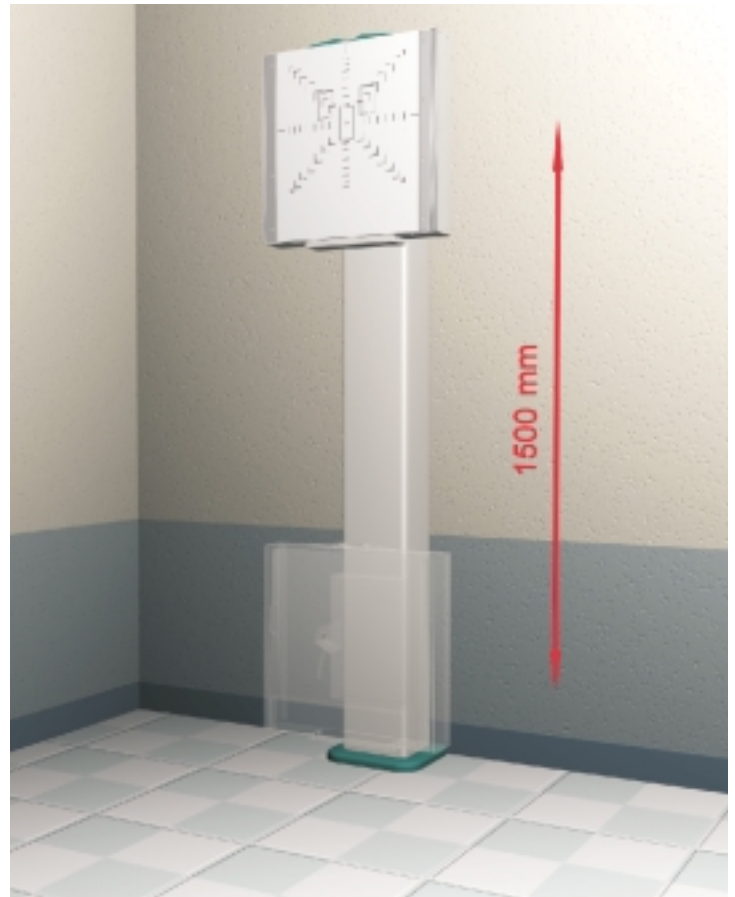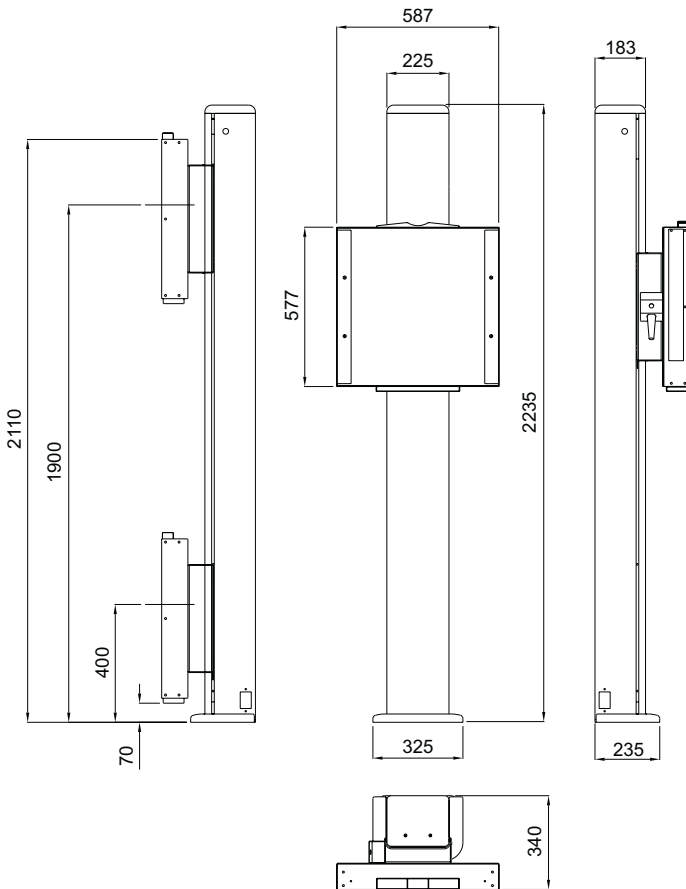

## Main Features

- Bucky vertical travel 150 cm
- Bucky center minimum distance from floor 40 cm
- Standard painting RAL 9002
- Weight 100 kg

## Accessories

- Arm support
- Compression band
- 35 x 43 cm external cassette holder
- 40 x 90 cm external cassette holder
- 40 x 120 cm external cassette holder

### SEDECAL

Pelaya, 9 - Poligono Industrial "Rio de Janeiro"  
28110-Algete (Madrid) - SPAIN  
Central Offices:  
Tel: (34) 91-628 0544 / 91-628 1592  
Fax: (34) 91-628 0574  
Email: marketing@sedecal.com

Factory:  
Tel: (34) 91-629 0721  
Fax: (34) 91-629 0783  
Email: manufacturing@sedecal.com

SEDECAL facilities



### SEDECAL USA Inc

2910 N Arlington Heights Rd.  
Arlington Heights, Illinois 60004  
USA  
Tel: 1 (847) 394 - 6960  
Fax: 1 (847) 394 - 6966  
Email: sales@sedecalusa.com